

All you need



MultiMode™ SuperPulse™ CO₂ Laser Platform

Your gateway to the future of aesthetic and women's health solutions



LUMENIS.COM/ACUPULSE

We don't just build technology, we create a world of solutions

AcuPulse™ is a complete ablative and resurfacing platform featuring state of the art technology for aesthetic resurfacing and women's health needs. Developed by Lumenis®, pioneer in energy-based technologies, it relieves a wide range of aesthetic and surgical issues as well as renewing the heart of a woman's being. AcuPulse™ restores not just health, but confidence allowing your patients to live life on their own terms, go forward and embrace everything it has to offer.

Live life to the fullest

Powerful. Versatile. Simple Platform.

SuperPulse[™] Technology

- Exceptional ablation/coagulation ratio
- Reach the desired treatment depth with a single pulse

Short Pulse Duration


- Enhanced patient comfort
- Minimizing thermal damage

11 Built-in Treatment Modes


- Maximum versatility with two scanners
- Full parameter control for advanced users

Intuitive User Interface

- Easy setup and quick treatments
- User presets for single-click treatment initiation

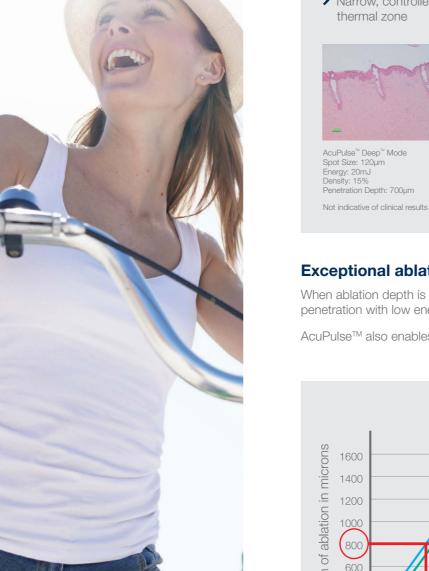


Unique SuperPulse[™] Technology And a Full Range of Advanced Continuous Wave Capabilities

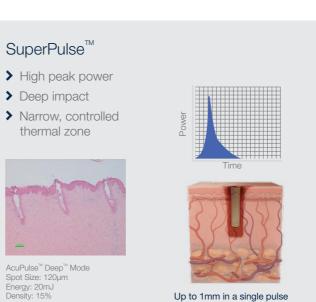


Exceptional ablation/coagulation ratio

When ablation depth is equal, more energy means more excess heat in the tissue. SuperPulse™ technology enables deep penetration with low energy. This reduces patient discomfort, downtime and thermal damage.

AcuPulse[™] also enables short pulses, further enhancing patient comfort.

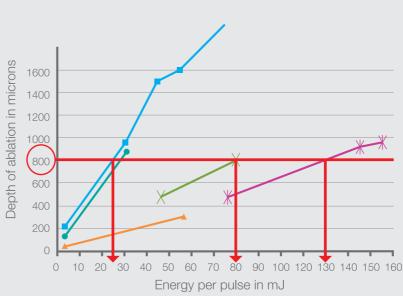

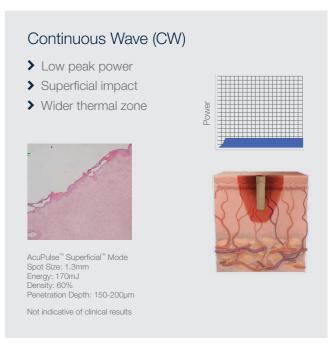





AcuPulse | 3

MultiModeTM: 11 Treatment Modes for Maximum Versatility

Aesthetic on AcuPulse[™]: A Trusted Path to Outstanding Outcomes



AcuScan120[™] fractional scanner

SurgiTouch™ advanced scanner



260mm, 200mm and 125mm handpieces for SurgiTouch[™]

Superficial[™]

For uneven pigmentation, fine lines, uneven texture, dyschromia, actinic keratosis and more.

Deep™

For deeper lines and wrinkles, acne scars, uneven texture, dyschromia and more.

Women's health

FeatherTouch[™]

For a variety of indication related to vulvar health

For gentle, single scan treatments.

SilkTouch™

For more vigorous, double scan treatments.

FineTouch™

For precise ablation of irregular pigmented lesions at various depths.

Combo[™]

Deep and Superficial combined together for uneven pigmentation, acne scars, wrinkles and more.

StretchTouch[™]

For skin furrows and textural irregularities.

ToeTouch™

For common foot ailments, such as plantar warts and ingrown toenails.

Paint[™]

For controlled, continuous freehand placement on small areas.

Women's health

For a variety of indication related to vulvar and vaginal health

Skin resurfacing

Courtesy of Rafael Nunes, MD



Skin resurfacing Courtesy of Girish Munavalli, MD







Blepharoplasty and lower lid resurfacing

Courtesy of Girish Munavalli, MD



Before

Before

After

Soft tissue (mole) removal Courtesy of Nariaki Miyata, MD



Before

After



125mm incisional handpiece

Bleph™

For precise, clean incisions with minimal bleeding.

"I have used the AcuPulse™ to perform over 500 CO₂ laser procedures. We have used it for skin conditions as varied as warts to syringomas to photo-damage. The Combo mode is especially helpful as the unique injury patterns allow for deeper fractional treatments (e.g. for acne scarring), while still allowing for lighter conventional superficial resurfacing for pigment and fine lines (all in one pass). Its SuperPulse[™] technology allows me to be more aggressive when I need to, while maintaining a wide patient safety margin, which ensures my patients are satisfied with the results of their treatments after a single session."

E. Victor Ross, MD, CA, US

"As a plastic surgeon, the AcuPulse $^{\rm TM}$ CO $_{\!\scriptscriptstyle 2}$ laser is perfect for me since I can offer effective treatment of scars. Many patients come to me because I treat with both CO₂ and ResurFX[™] lasers soon after surgical procedures. They are confident they won't have scars, and they don't."

Francesca de Angelis, MD

Women's health on AcuPulseTM: 360° of women's health solution

The women's health on AcuPulse[™] offers a range of minimally invasive solutions for women coping with gynecologic conditions that affect their well-being, confidence and intimacy.

Gynecological treatment

When operating on delicate anatomy, the aim is, first and foremost, to avoid thermal damage. That's why the AcuPulse™ CO, laser technology is a treatment of choice, suitable for use with a variety of approaches and applications.

Vaginal approach

Colposcopy approach

Cervical Intraepithelial Neoplasia (CIN)

- ✓ Vaginal abscess
- Genital warts
- Labial reduction
- Lichen sclerosis
- Leukoplakia Condyloma
- Benign & malignant tumors

Laparoscopy approach

- Endometriosis
- Adhesions
- Myomas
- Cysts
- Hysterectomy

Your journey with Lumenis has just begun.

We want your business to grow with Lumenis That's why we have put together a support package that begins with the first that you have all the training you need and receive the most up to date clinical materials.



Lumenis Practice Locator

Thousands of patients search for treatment providers every month. Help them find your clinic by registering your practice on the Lumenis Practice Locator. We run campaigns through influencers, social media and viral videos to promote this directory.

Lumenis online

Consumers have many questions about aesthetic. Our patientfacing websites will give them the re-assurance and confidence they need to move forward with the treatments.

000

≣



Your Most Valuable Resource

Your dedicated representative is an email, text, or phone call away. We want you to feel part of the Lumenis family, so please feel free to contact us any time. We'll always do our best to get back to you within 24 hours.

Brand Your Clinic

As a Lumenis client you will have access to the Lumenis Physician Portal. Here, you can find materials such as brochures and posters. You will also have access to a wealth of clinical and product information.

Team Training

We install the system and train your staff, so that you are able to assess skin types and conditions, creating the perfect treatment plan for each individual patient.

Our expert clinical training means they will handle consultations and treatments with confidence.



Extensive Service Network

We are committed to providing high quality, comprehensive service for medical equipment repair and maintenance. Lumenis service provides you with peace of mind, so you can focus on patient care, rather than on the equipment.

Specifications

Platform	AcuPulse [™] 40AES-F, 40AES-R, 40AES-A
Technical Specifications	
Wavelength	10,600 nm
Laser Type	SuperPulse [™] and CW
Power to Tissue	40 W
Peak Power	170 W
Spot Size	min. 0.12 mm / max. 1.3 mm (AcuScan120™)
Scan Size	Up to 15 x 15 mm
Density per Scan	fractional: 3-60%, full ablative: >100%
Full Resurfacing Capabilities	yes
Fractional Capabilities	yes
Incisional Capabilities	yes
Excisional Capabilities	yes
Aiming Beam	5 mW red diode laser, 635 nm, adjustable intensity
Beam Delivery	Lightweight, 7-joint, spring-balanced arm w/120 cm (47.2 in) reach, 360° rotation
Memory Setting	min. 100
Cooling	Self-contained, closed cycle
Warranty and Coverage	1 year standard
Recommended Passes/Treatment	1
Platform Specifications	
Electrical	100–240 VAC, 10 A max., 50/60 Hz
Weight	49 kg (108 lbs)
Dimensions (WxDxH)	37 cm x 40 cm x 119 cm (14.6 x 15.8 x 47 in)

Warnings and risks

CO₂ lasers are intended solely for use by physicians trained in the use of the Carbon Dioxide laser (10.6 µm) wavelength. Incorrect treatment settings or misuse of the technology can present risk of serious injury to patient and operating personnel. Risks that may be associated with any CO₂ laser procedure may include change of pigmentation, infection, erythema, skin induration or scarring. Read and understand the CO₂ systems and accessories operator manuals for a complete list of intended use, contraindications and risks.

The use of Lumenis® CO2 laser is contraindicated where a patient has taken Accutane (Isotretinoin) within the past 6-12 months, has a history of keloid formation and demonstrate excessive or unusually prolonged erythema.

© 2018 All Rights Reserved. The Lumenis Group of Companies. Lumenis, its logo, AcuPulse, UltraPulse, SuperPulse, AcuScan120, SurgiTouch, FeatherTouch, SilkTouch, FineTouch, ToeTouch, StretchTouch and Combo are trademarks or registered trademarks of the Lumenis Group of Companies. Specifications are subject to change without notice. Lumenis is ISO 13485:2003 certified. AcuPulse is CE approved. AcuPulse is FDA cleared. PB-1135190 Rev G.



Lumenis Ltd. Yokneam Industrial Park 6 Hakidma Street P.O.B. 240 Yokneam 2069204, Israel Tel: +972-4-959-9000

Lumenis (Germany) GmbH Heinrich-Hertz-Str 3 D-63303 Dreieich-Dreieichenhain GERMANY Tel: +49 (0) 6103 8335 0

EMEA Dreieich Dreieichenhain, Germany T +49 6103 8335 0

F +49 6103 8335 300 Roma, Italy T +39 06 90 75 230 F +39 06 90 75 269

Hertfordshire, UK T +44 20 8736 4110 F +44 20 8736 4119

Tokyo, Japan T +81 3 6743 8300 F +81 3 6743 8301

JAPAN

AMERICAS

San Jose, CA, USA T +1 408 764 3000 +1 877 586 3647 F +1 408 764 3999

ASIA PACIFIC Beijing, China

T + 86 10 5737 6677 F +86 10 5737 6767

Gurgaon, India T + 91 124 422 07 95 Kowloon, Hong Kong T + 852 217 428 00 F + 852 272 251 51

CE 0483

PB-1135190 Rev H

© LUMENIS LTD 2018



f 🖌 📾 in

LUMENIS.COM/ACUPULSE

For more information: www.lumenis.com